

leggera
FIRST STEPS GUIDE


STEP 1
UNVEIL YOUR SCALE
p.04

STEP 2
PAIR YOUR LEGGERA
p.08

STEP 3
PERFORM YOUR FIRST MEASUREMENT
p.08

STEP 1


UNVEIL YOUR SCALE.




A. Unpack your **Leggera** and check all the contents.



Smart Bluetooth
Leggera Scale

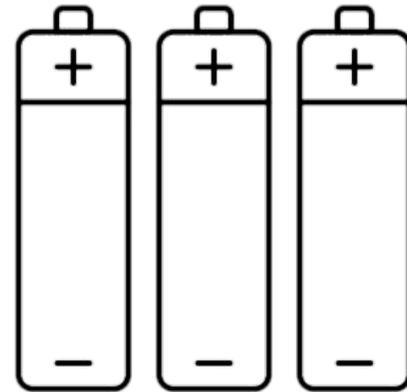



3x1.5V AAA
Batteries

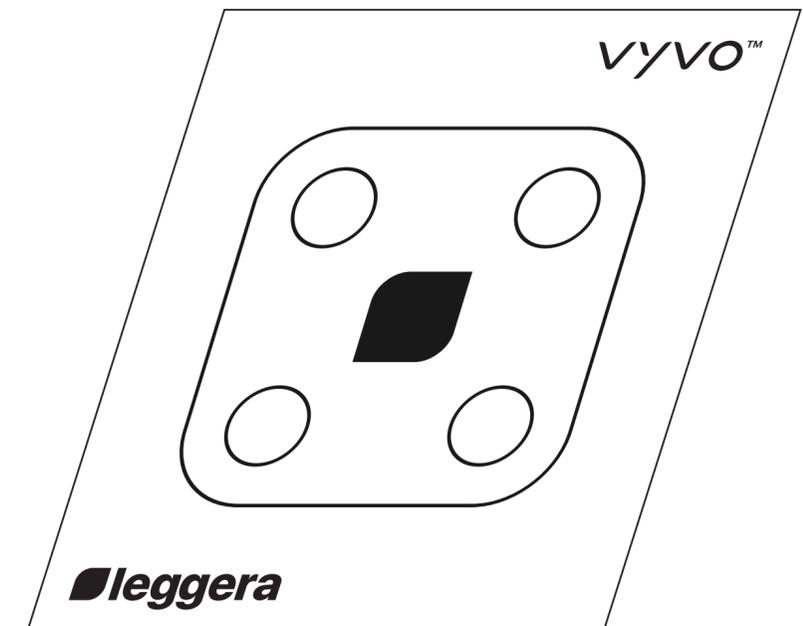

Quick Start
Manual



B. Put the batteries inside the Scale.

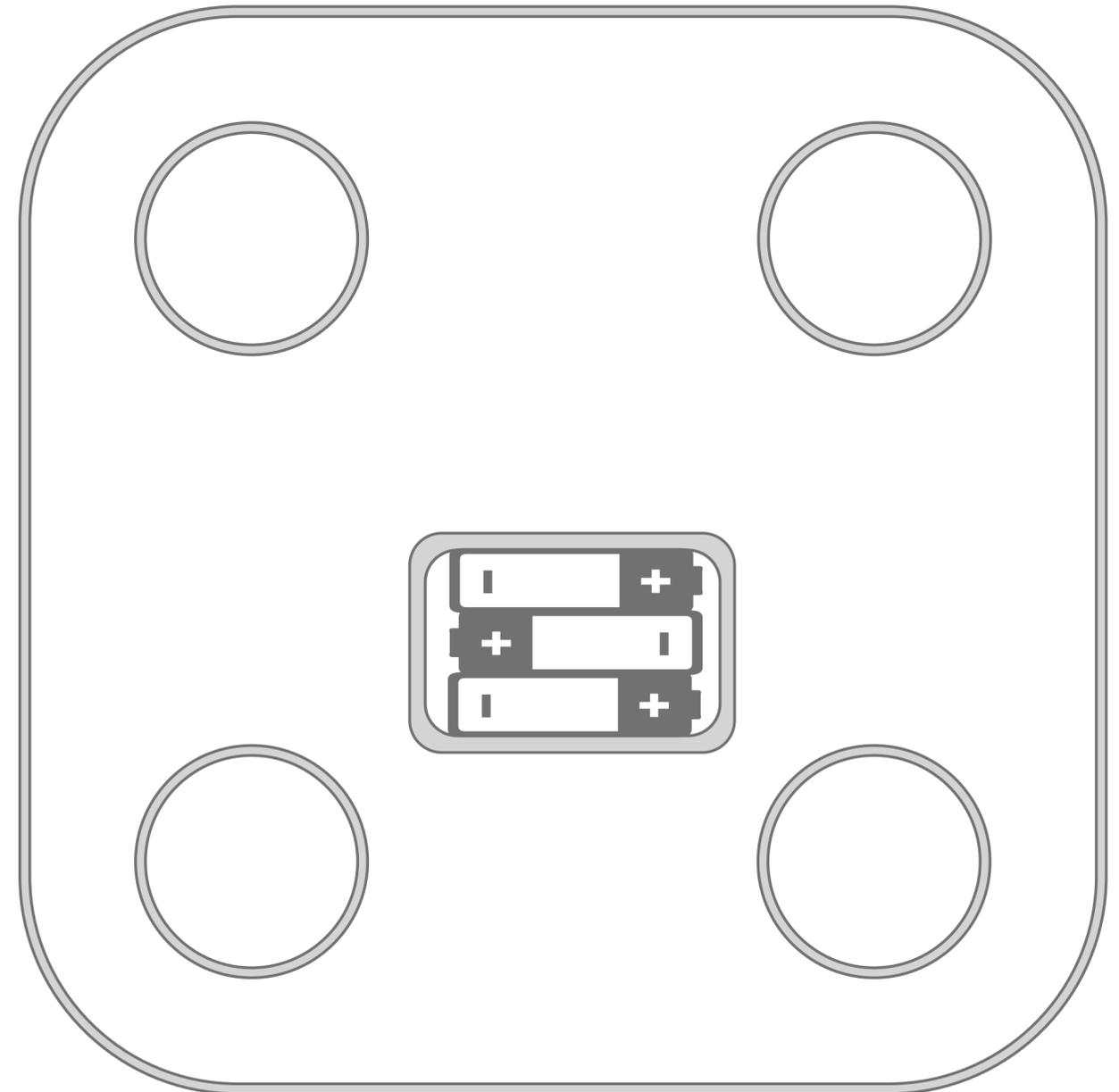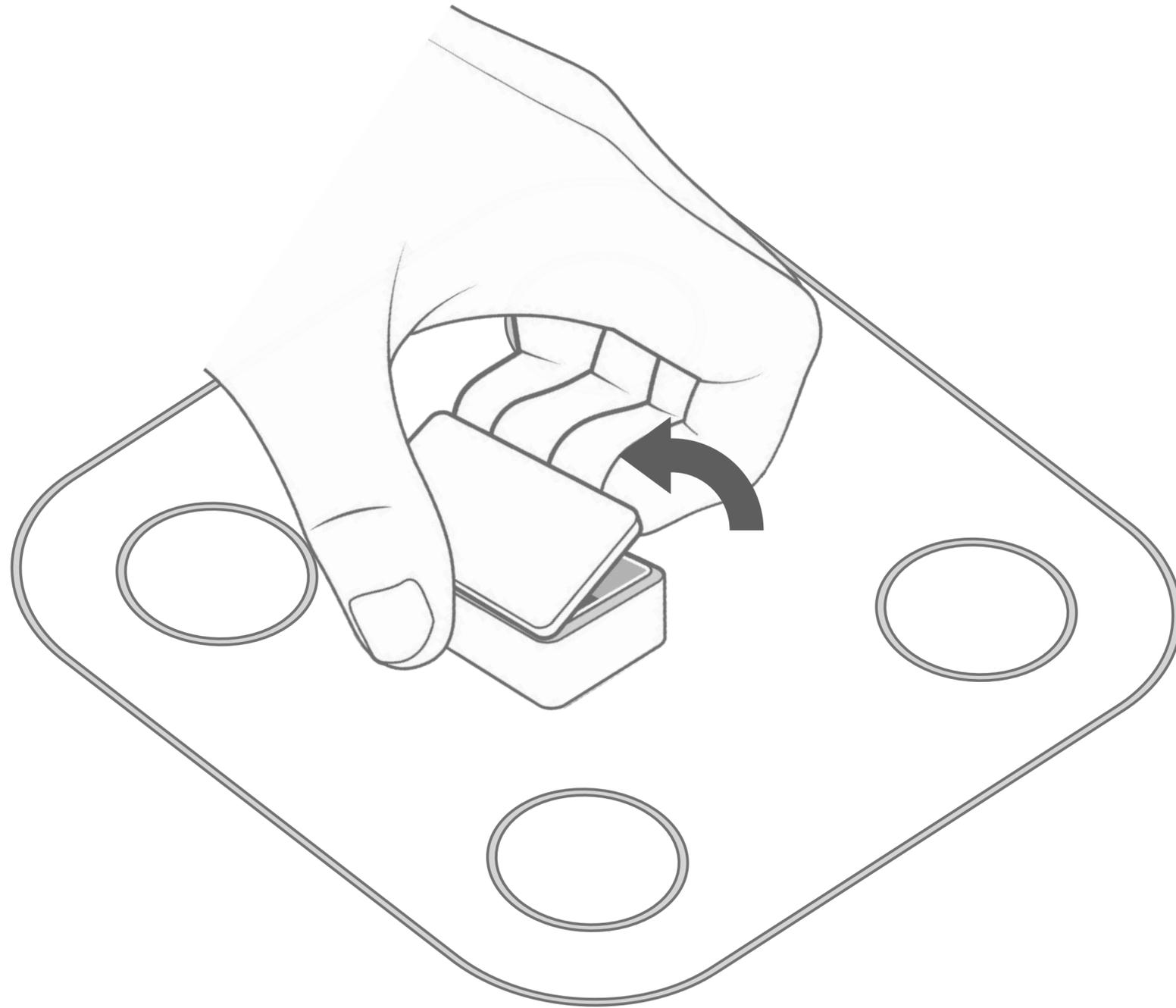

C. Place Leggera on the floor.





 **leggera**

STEP 2

PAIR YOUR LEGGERA.



A. Open the VYVO Smart App and access the Device section - Tap on "+".

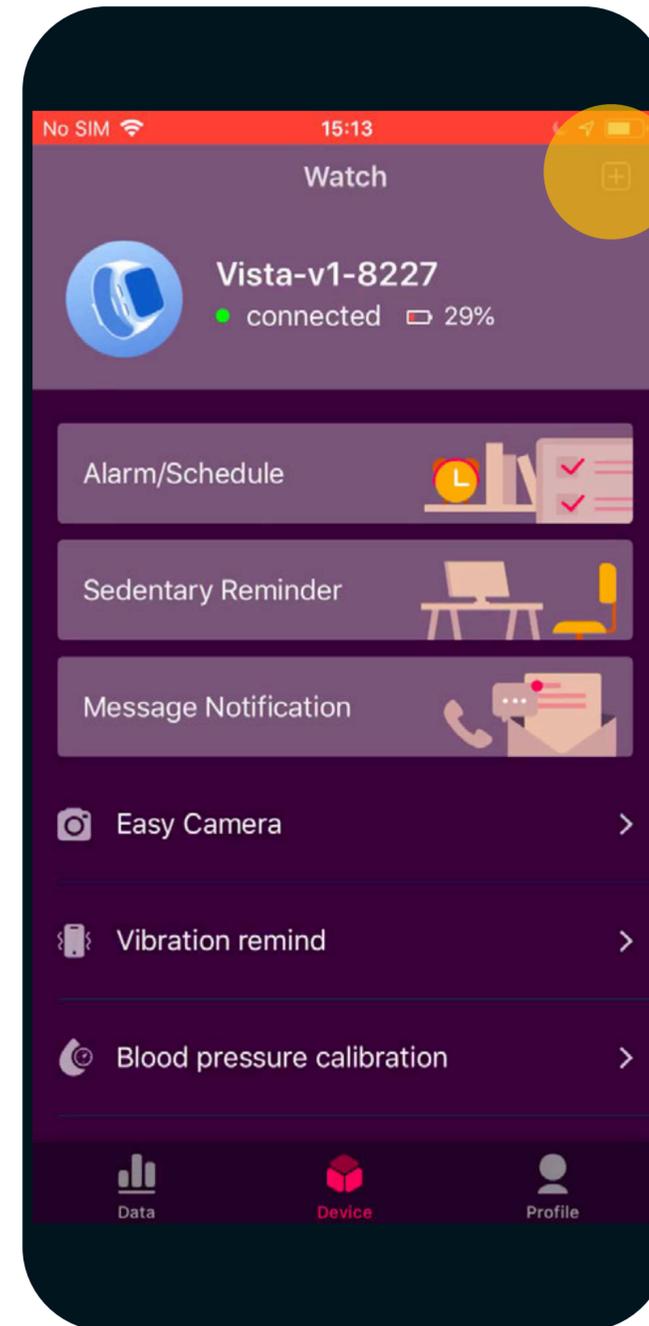
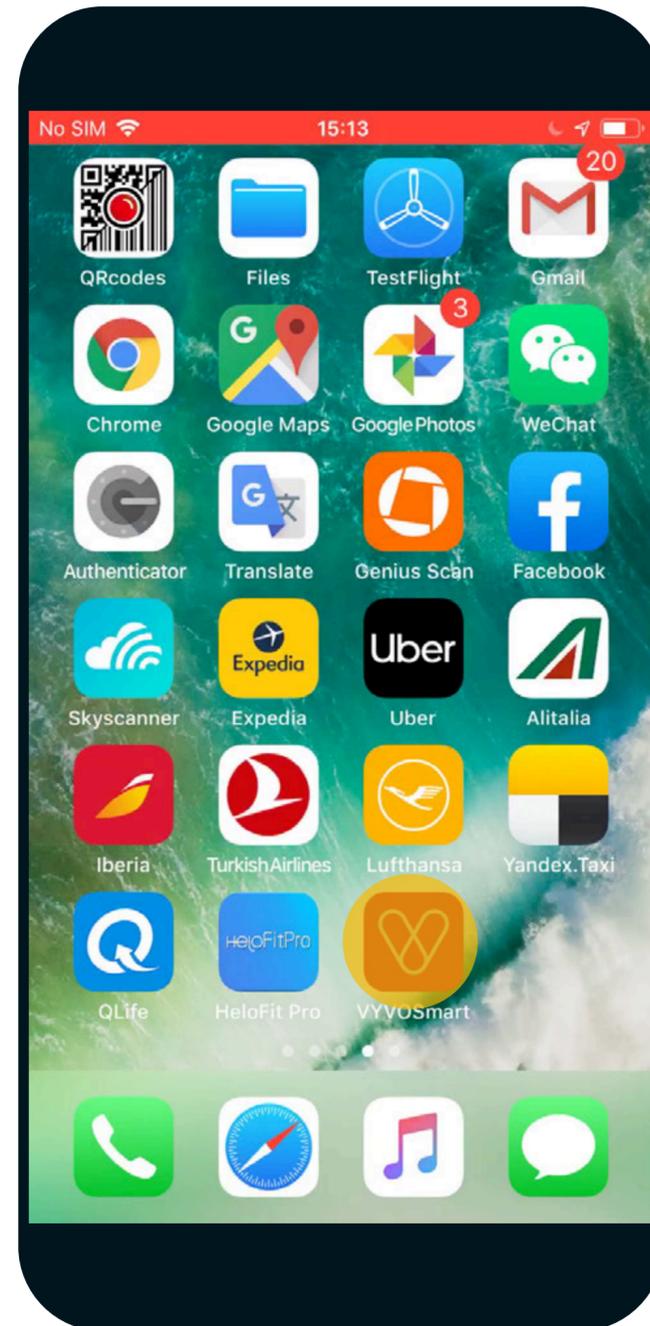


B. Remove your shoes and Step on Leggera to turn it on.





 **leggera**



C. In the VYVO Smart App tap on Health Scale and wait for pairing...

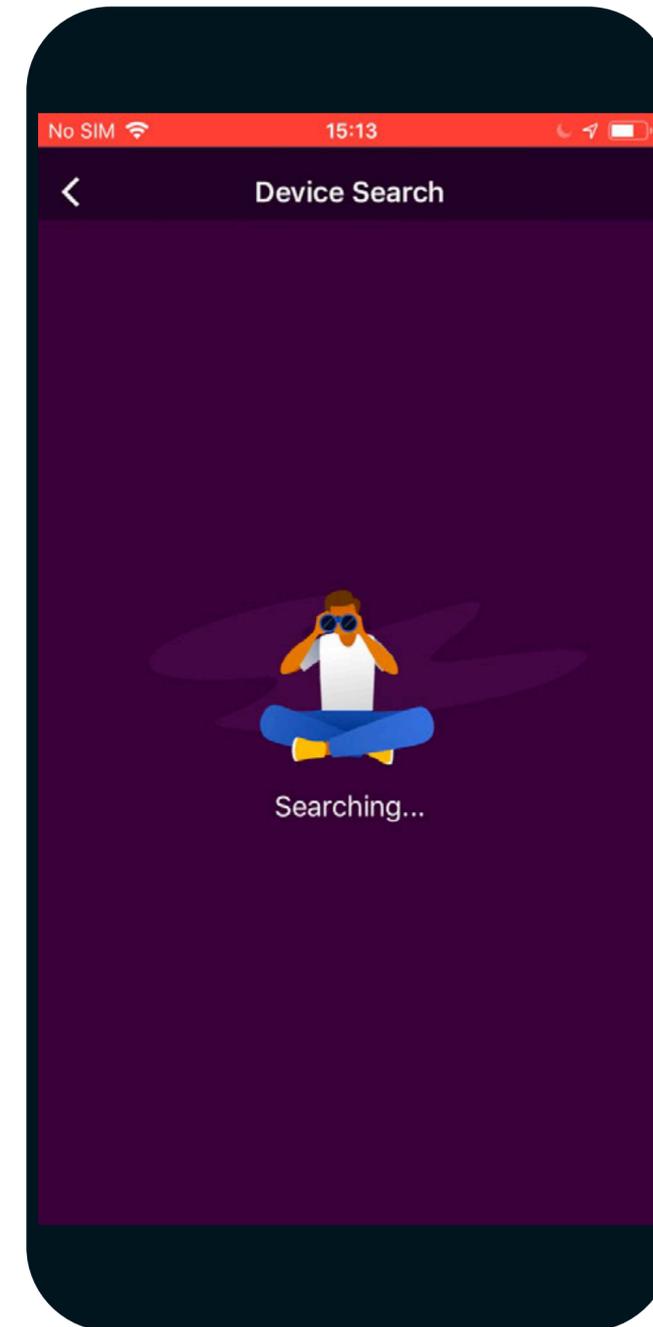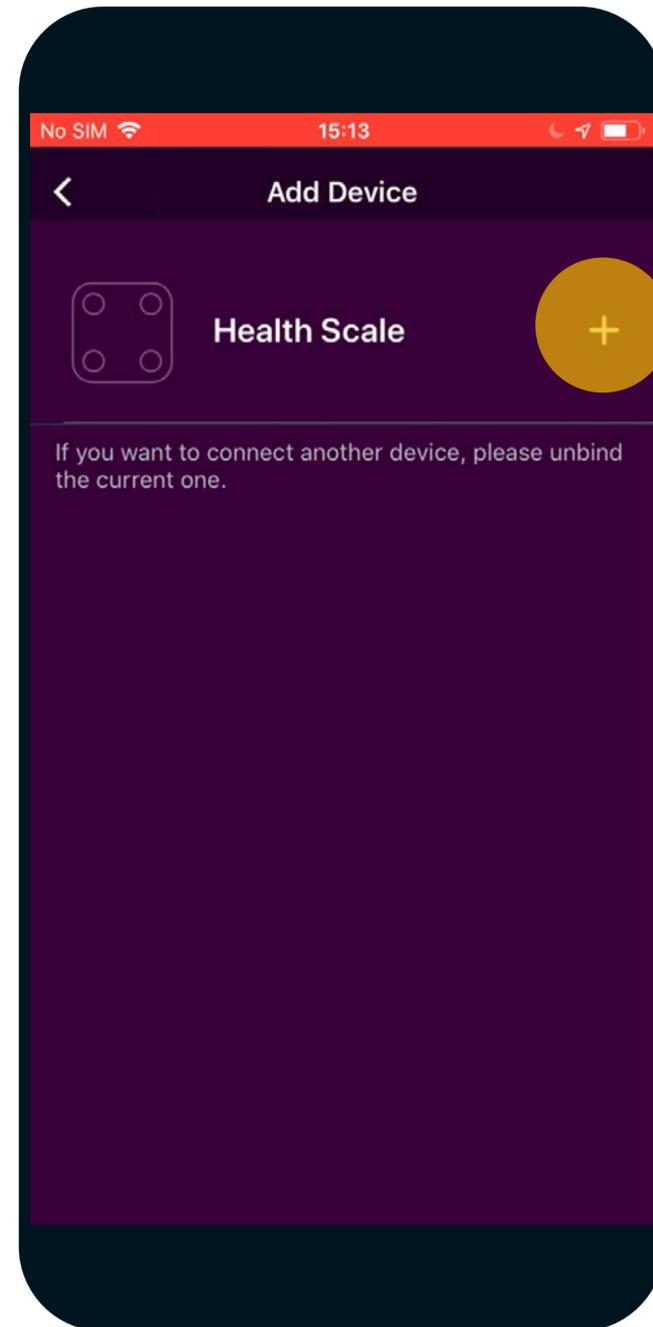




D. Click to add your Leggera and wait for the connection to complete.

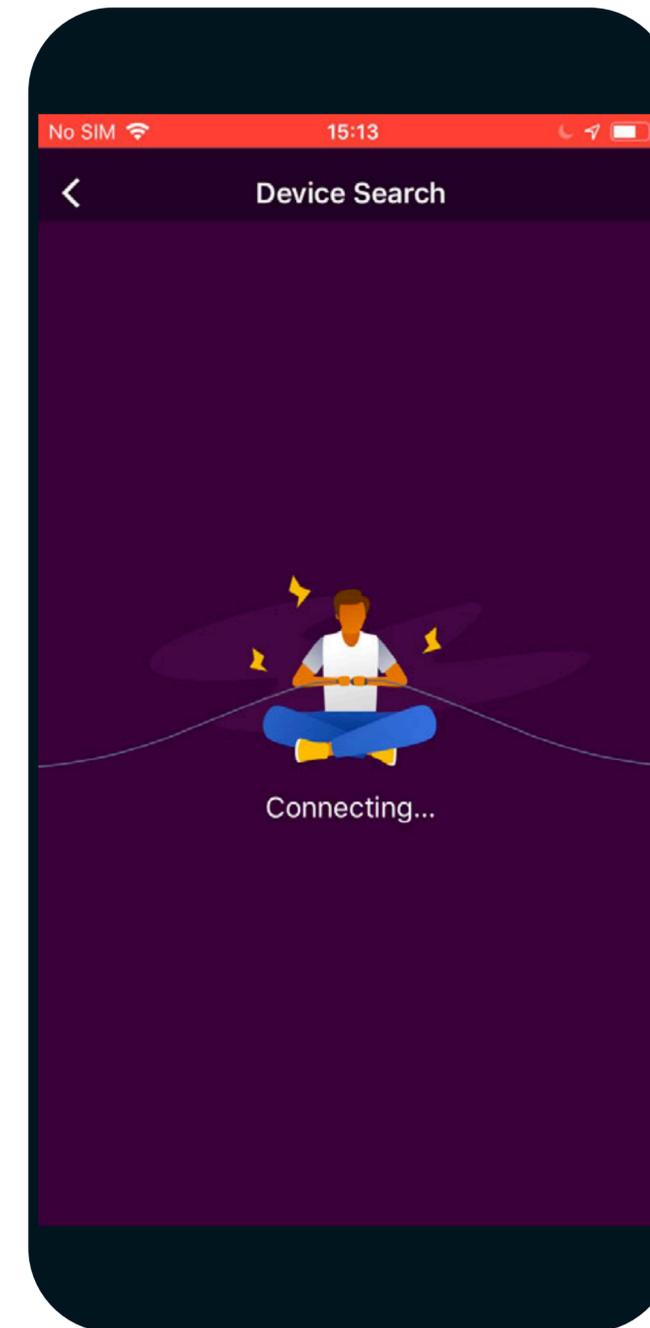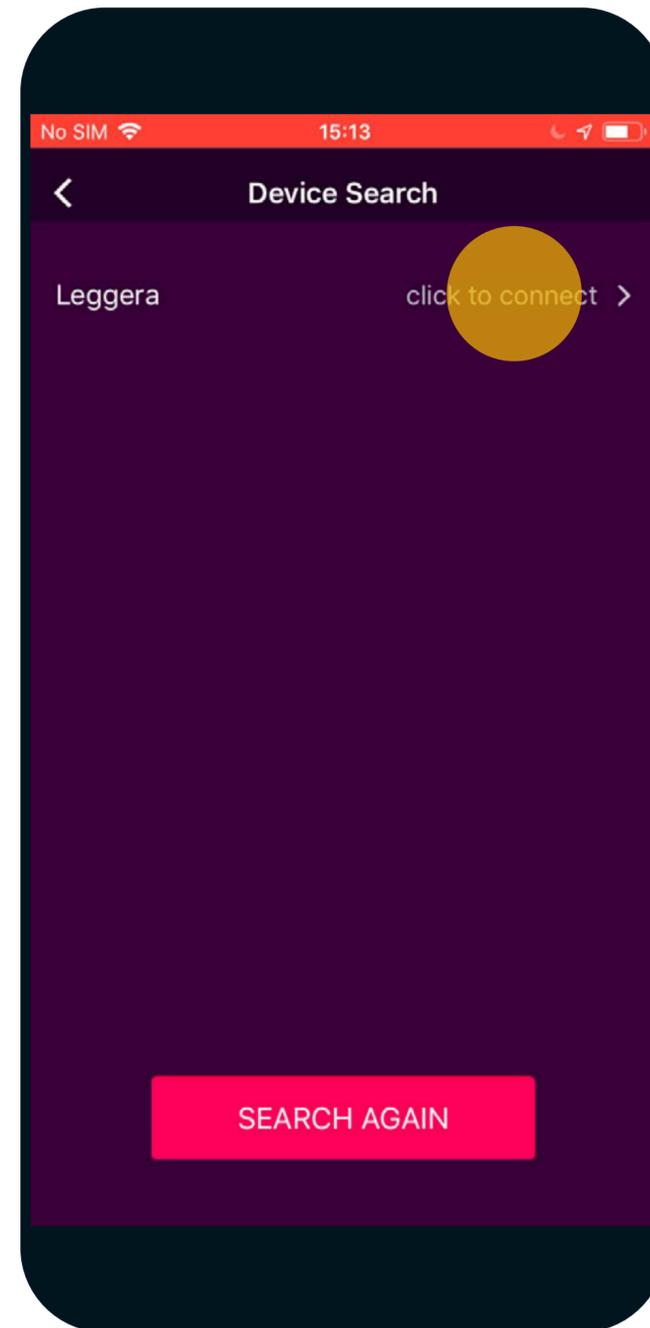




E. Set your preferred Leggera units of measurement.

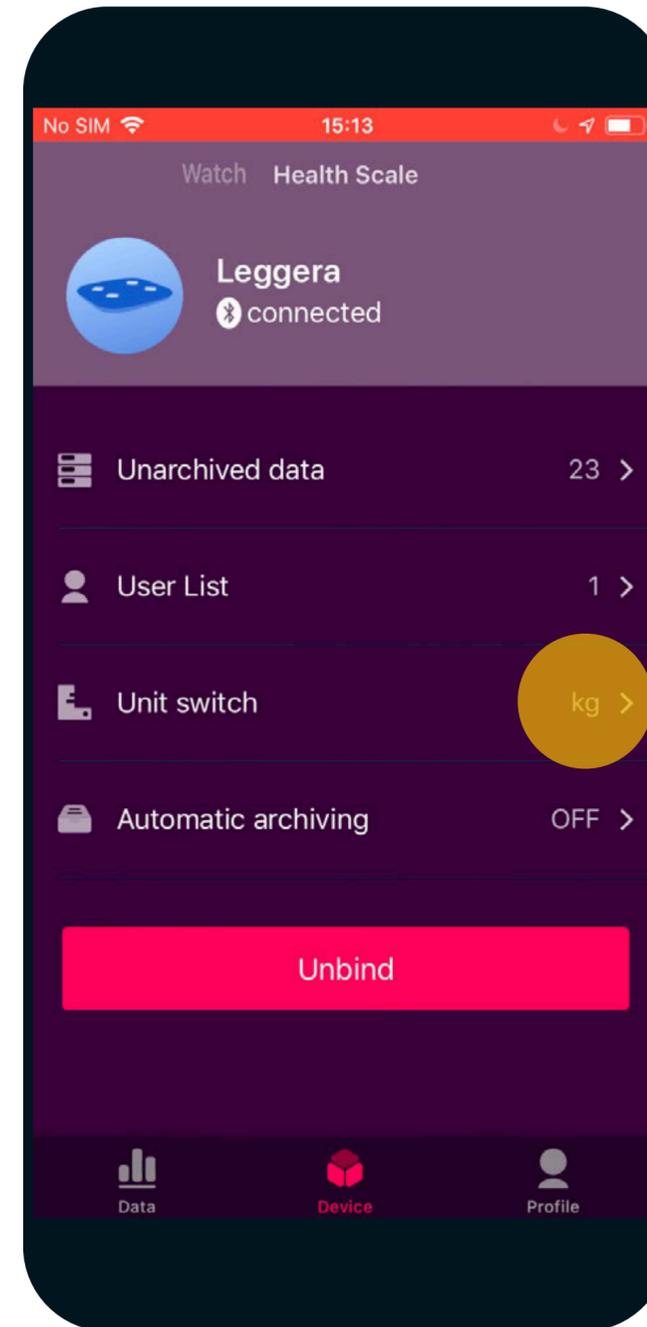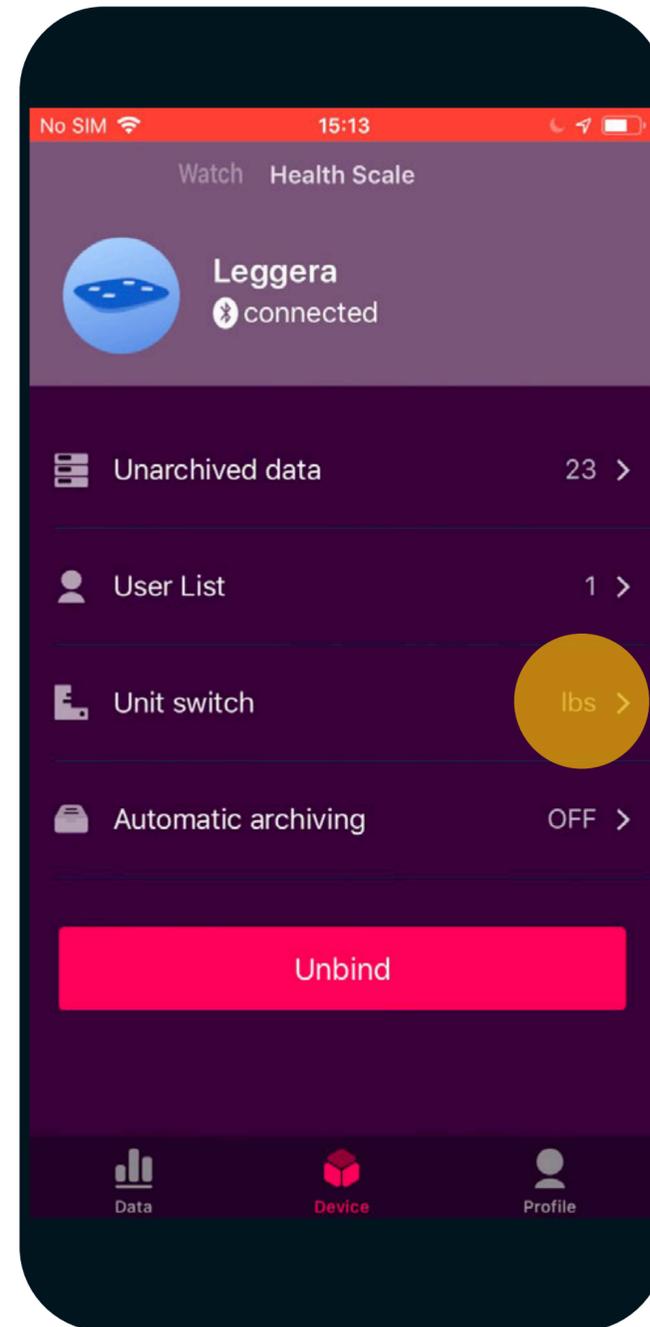




F. Set your preferred App units of measurement within "Profile".

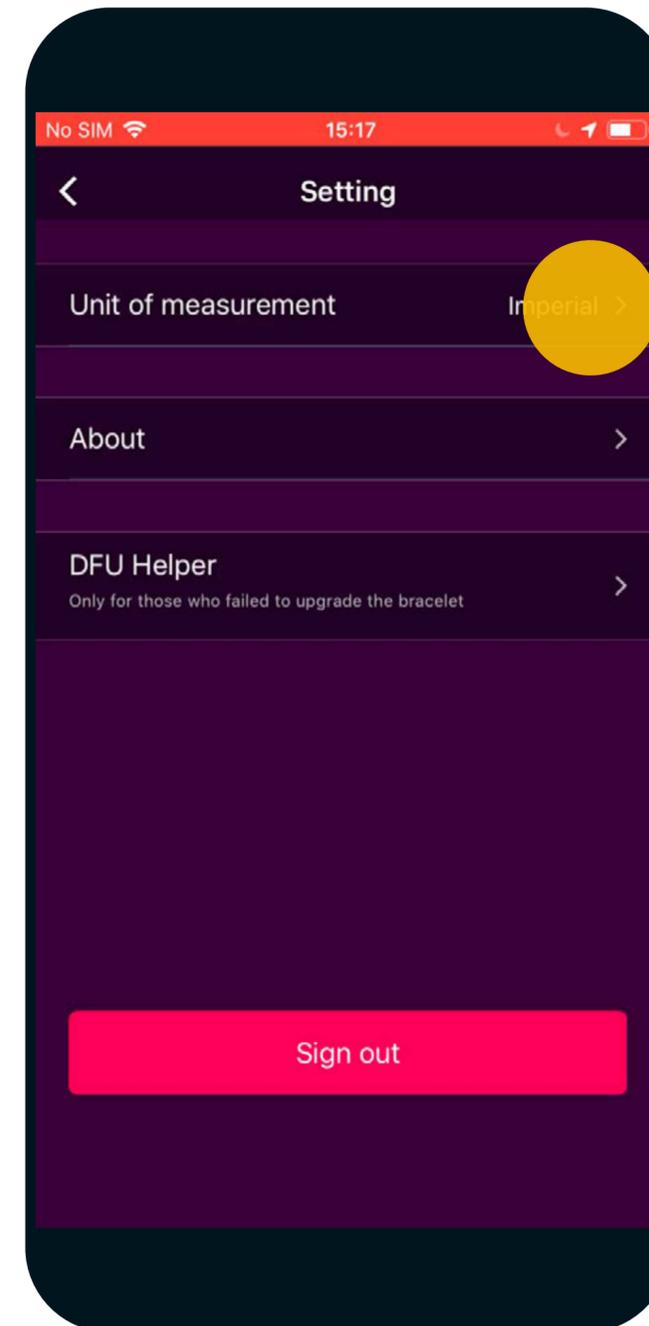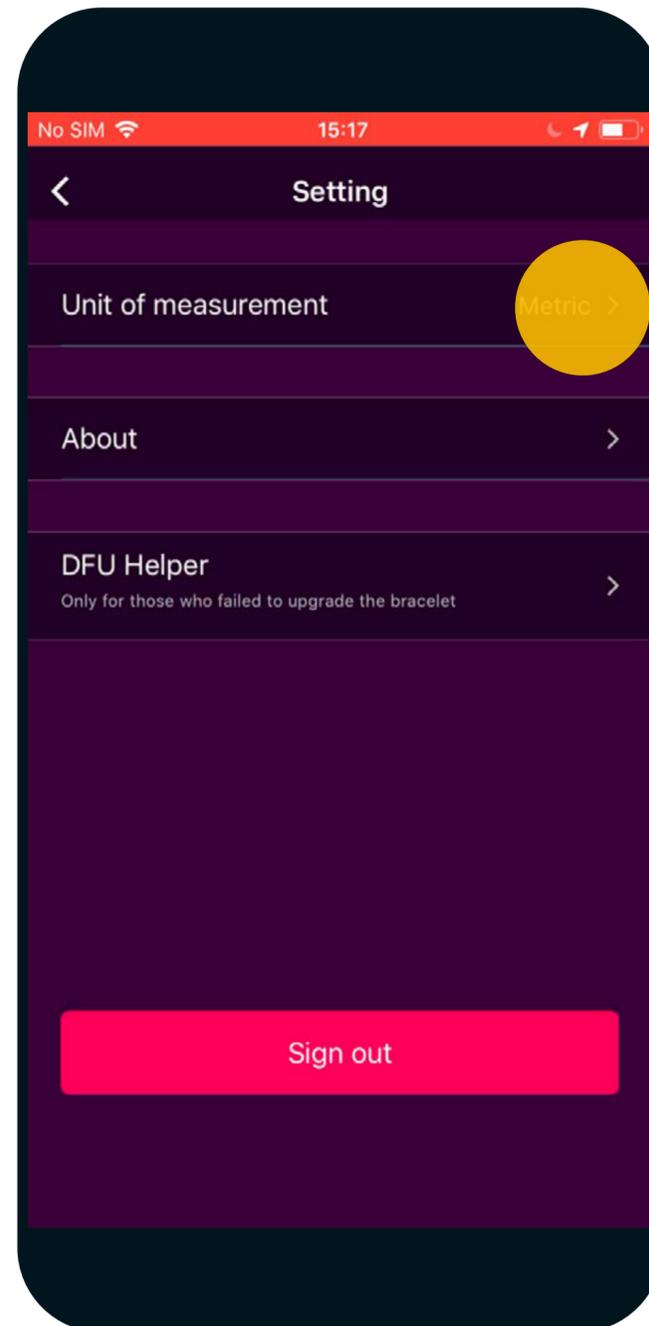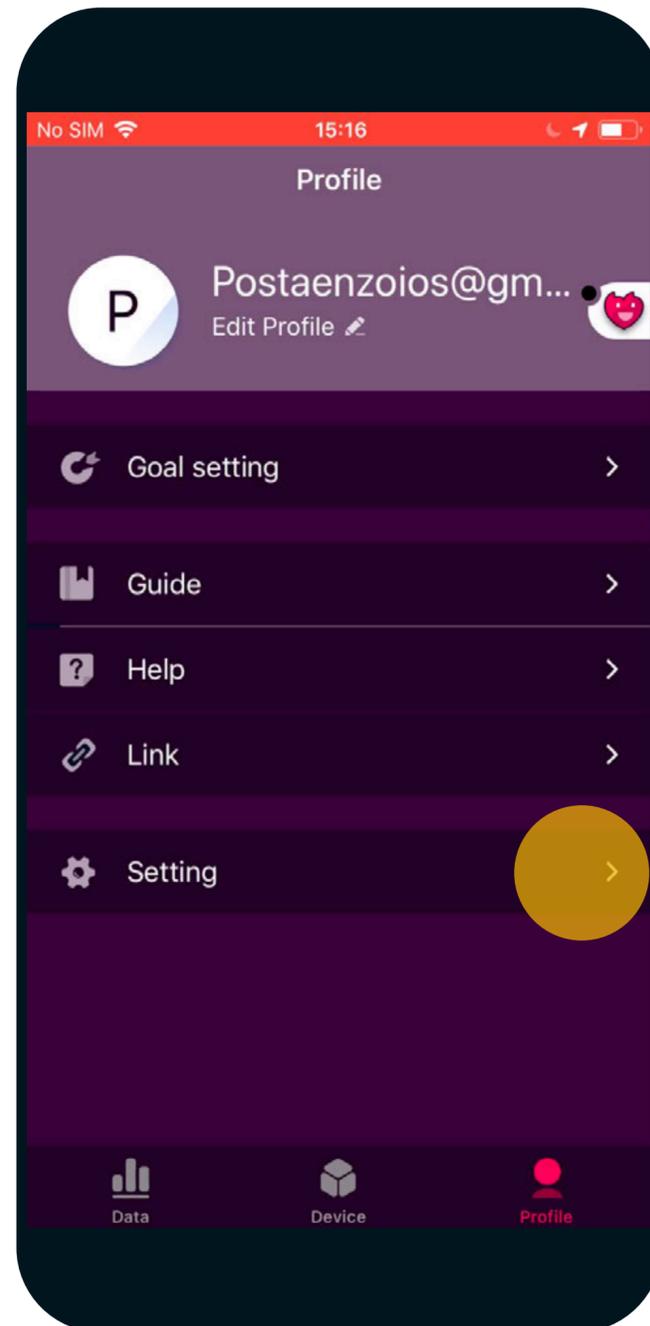

G. Add a new user profile.

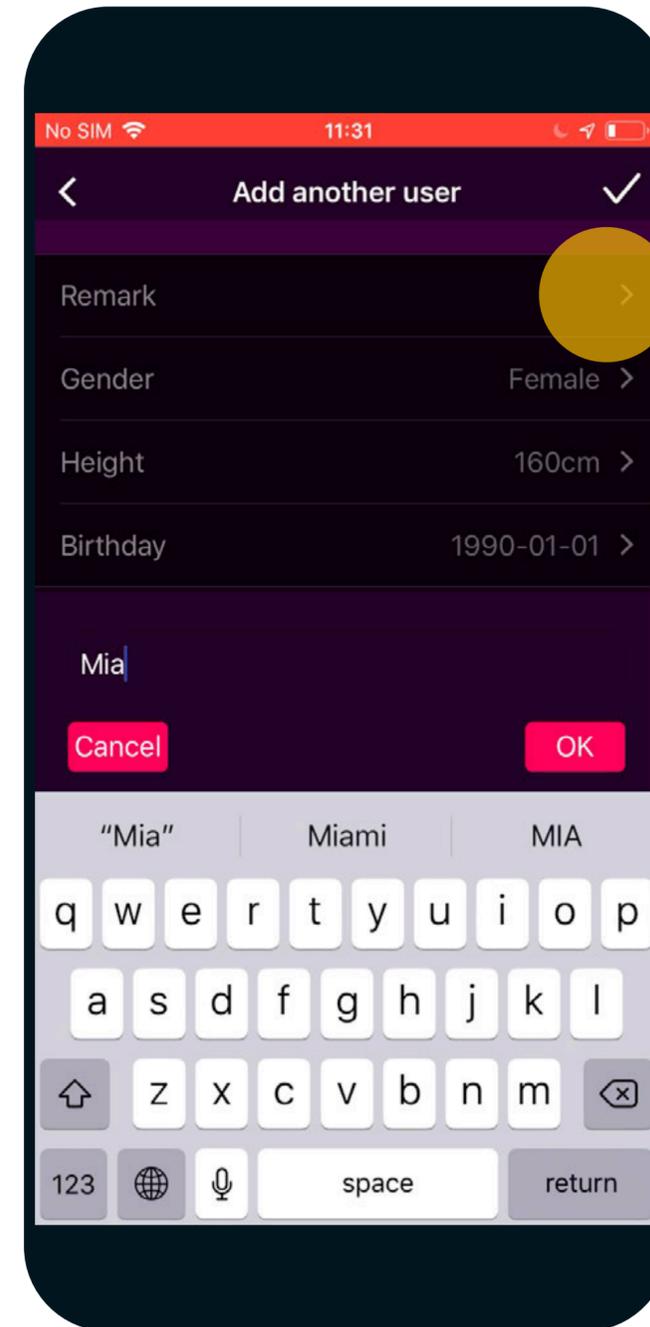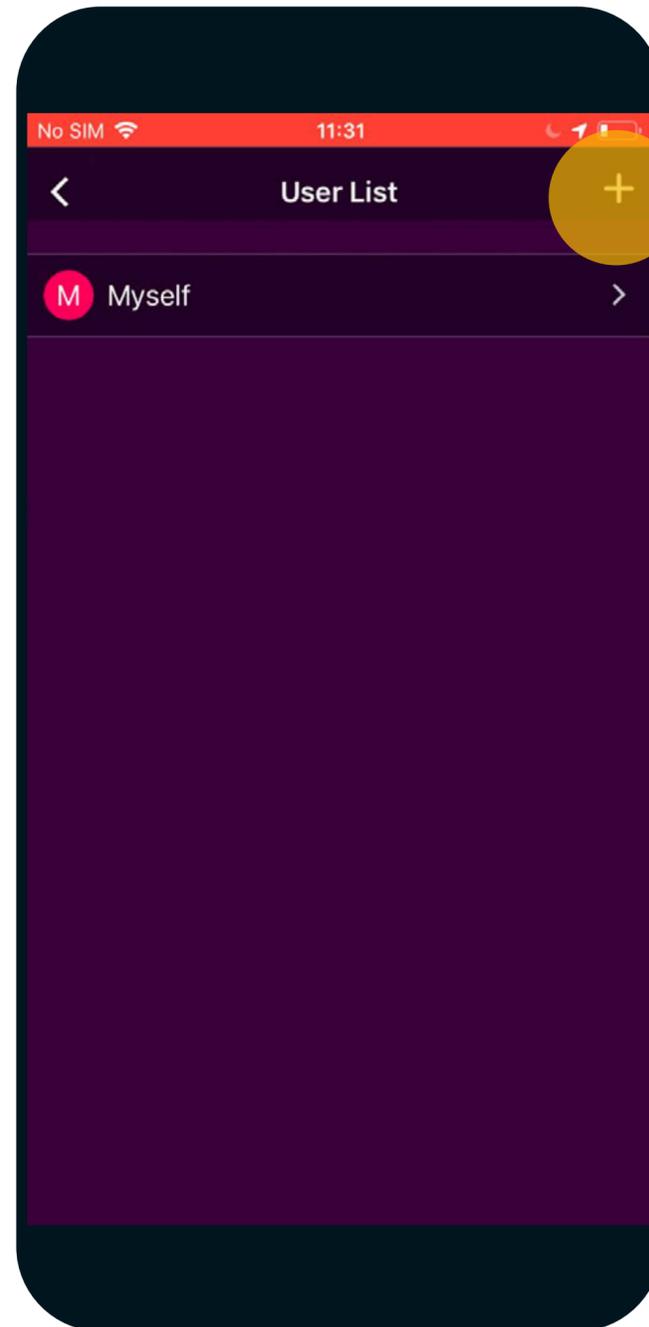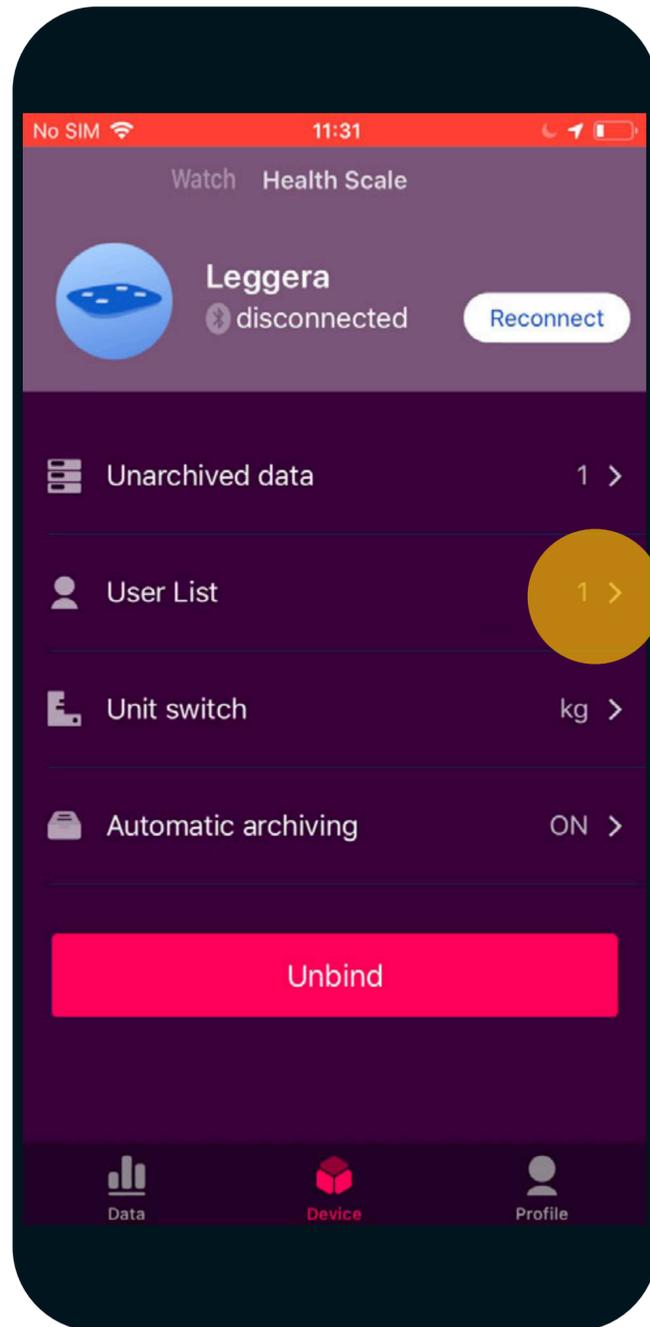

G. Add a new user profile.

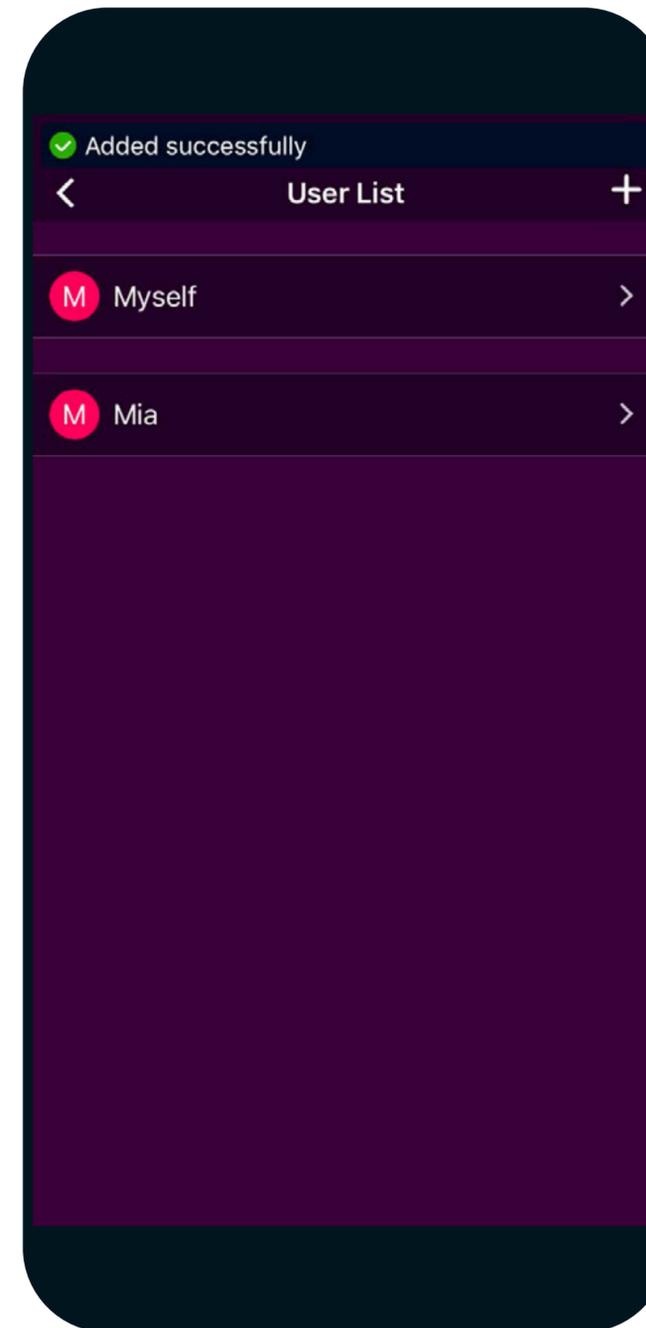
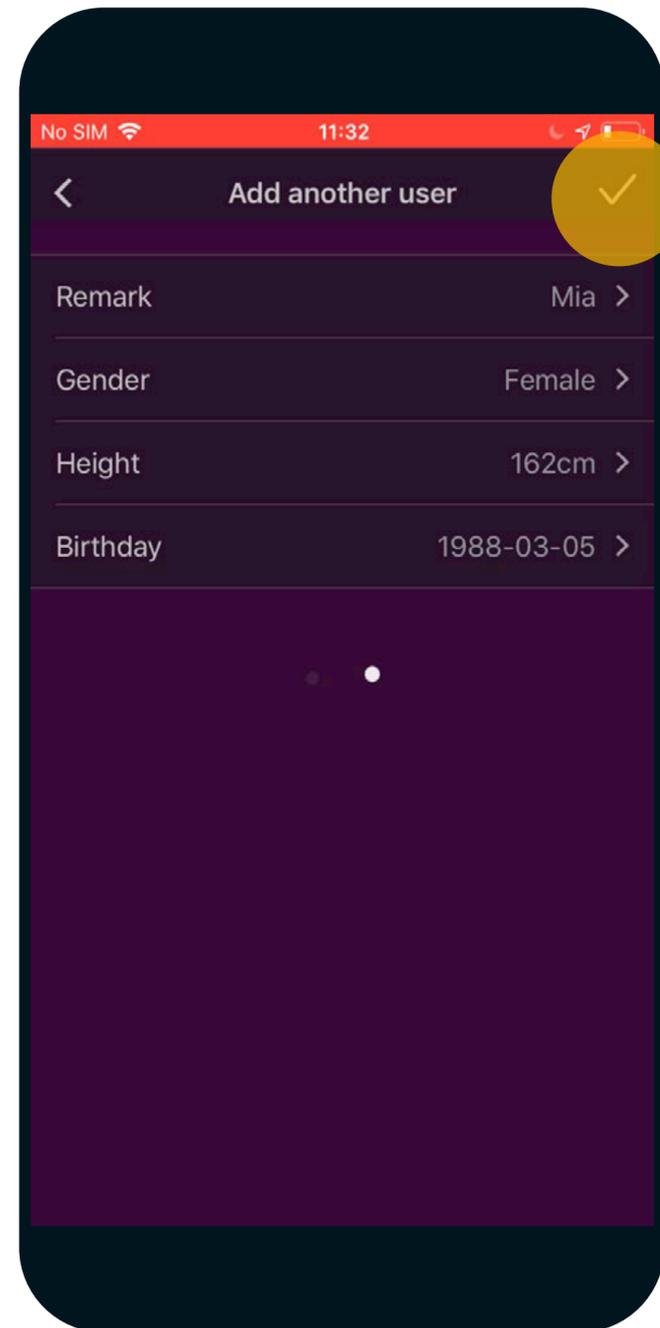
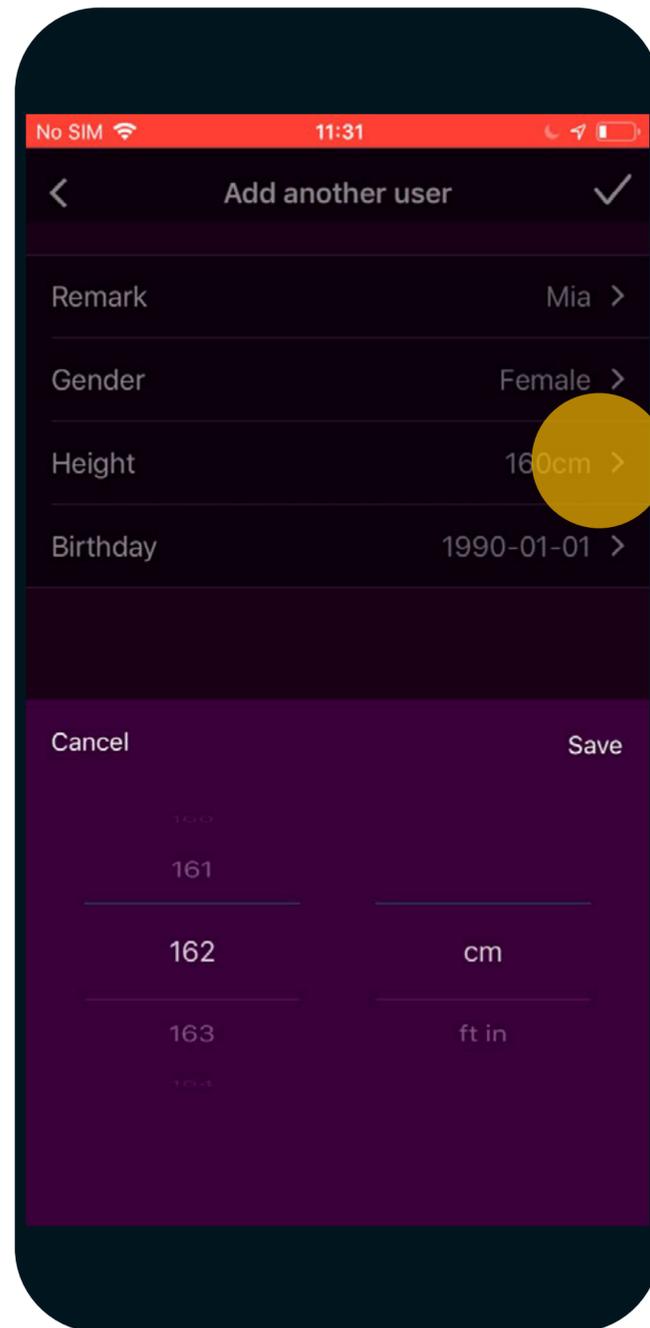


STEP 3

PERFORM YOUR
FIRST MEASUREMENT.



A. Open the VYVO Smart App and access the Device section Health Scale.

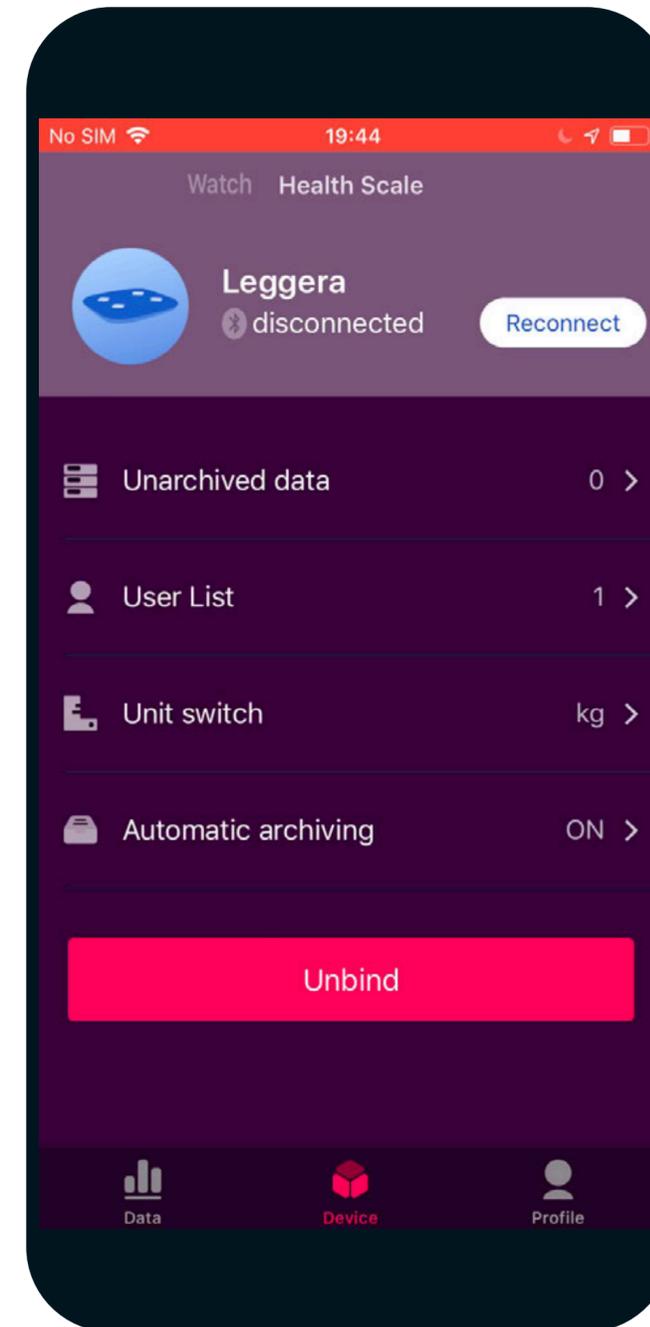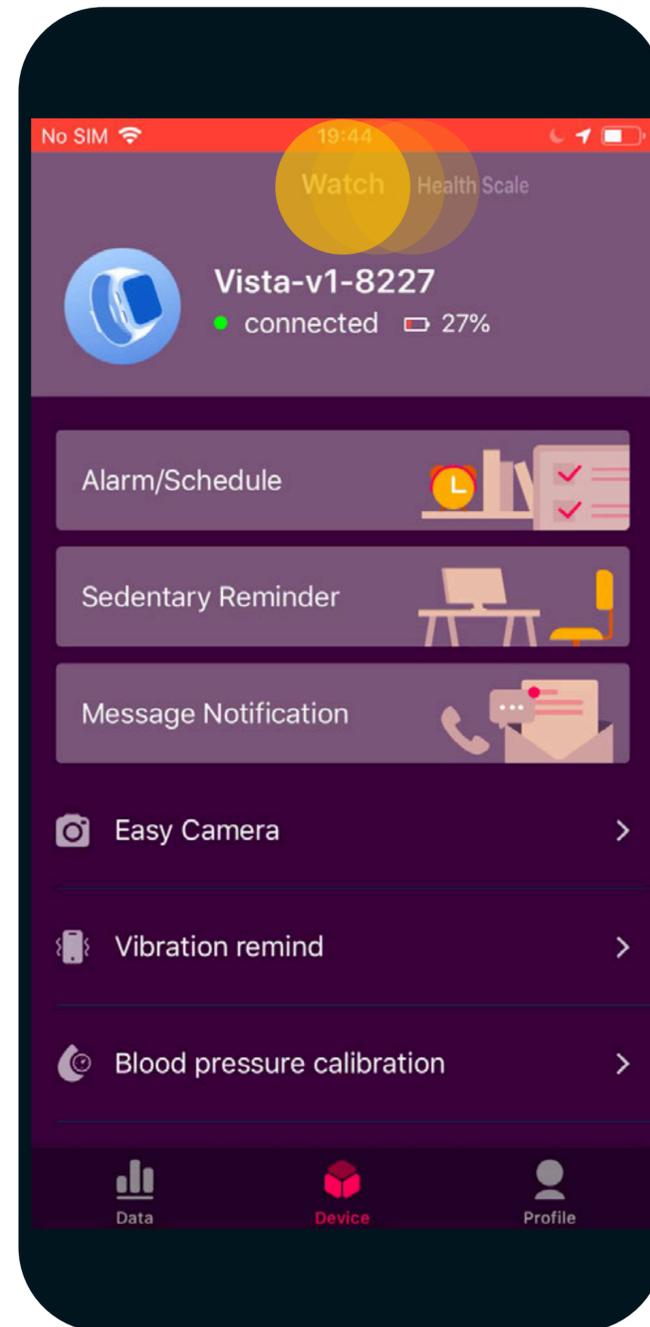
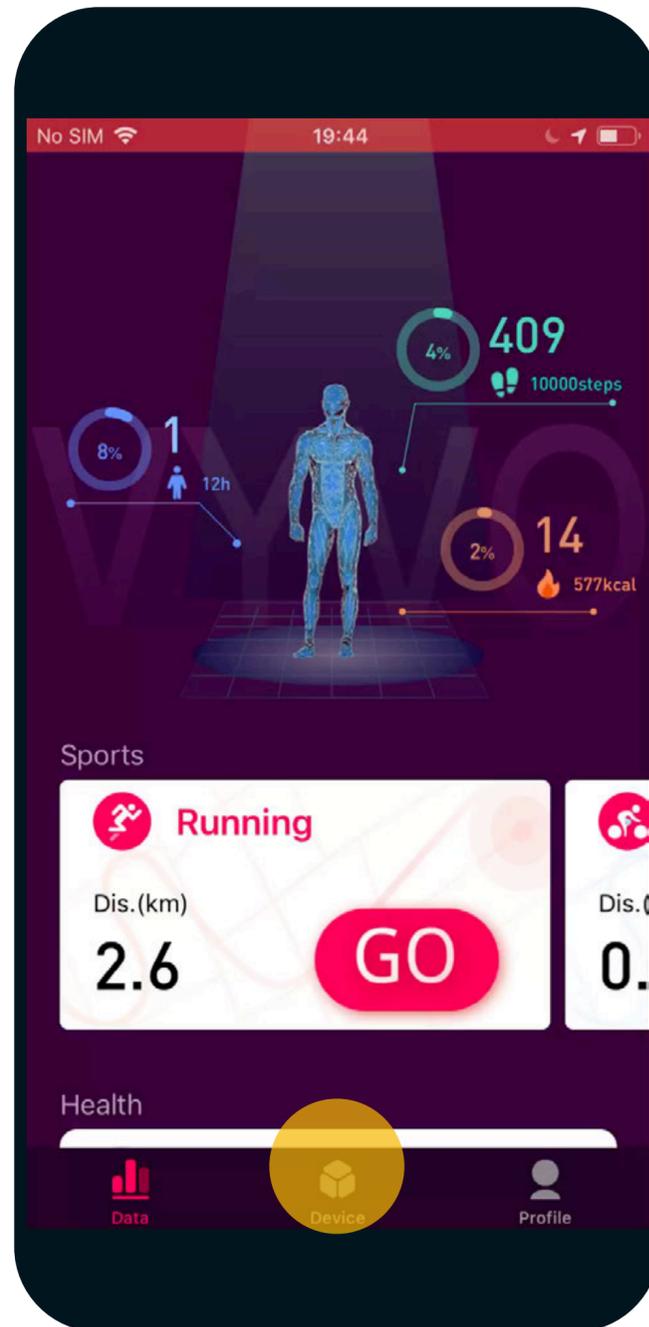


Swipe the screen on the left if your App is paired with a VYVO Watch, to access the Health Scale section.


B. Remove your shoes and Step on Leggera to start the measurement.




C. When the measure is completed, the unarchived data section will be updated and the count will increase.

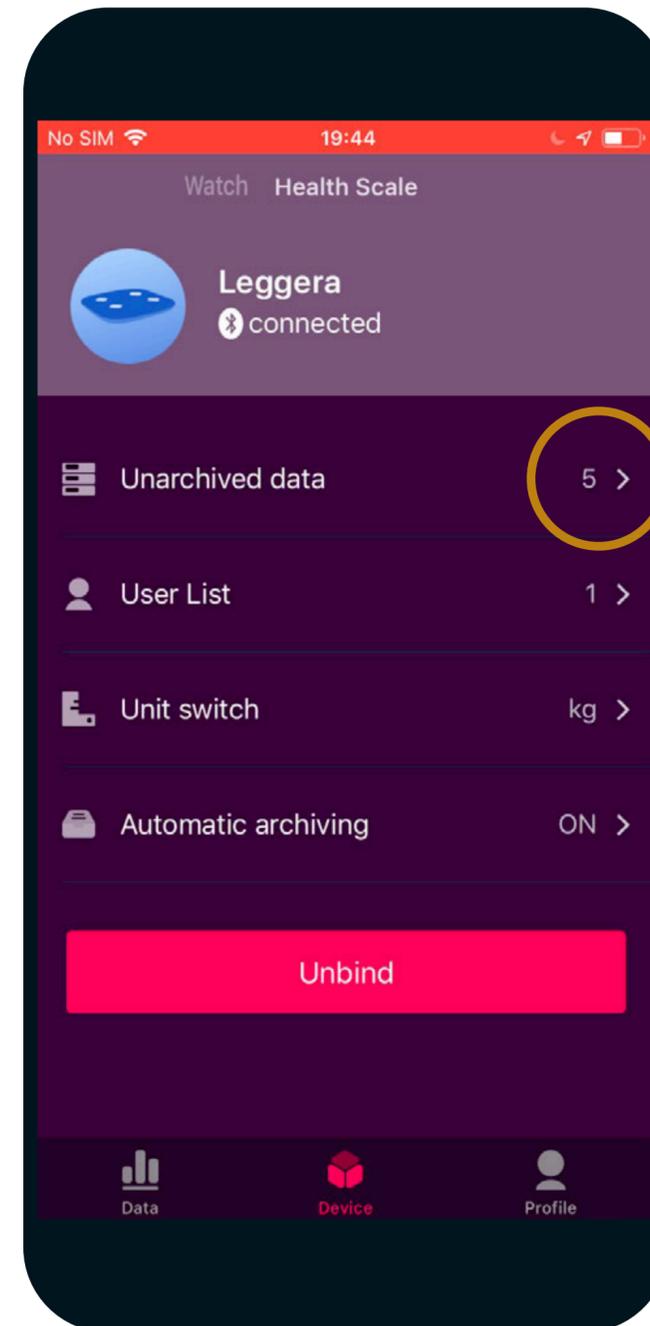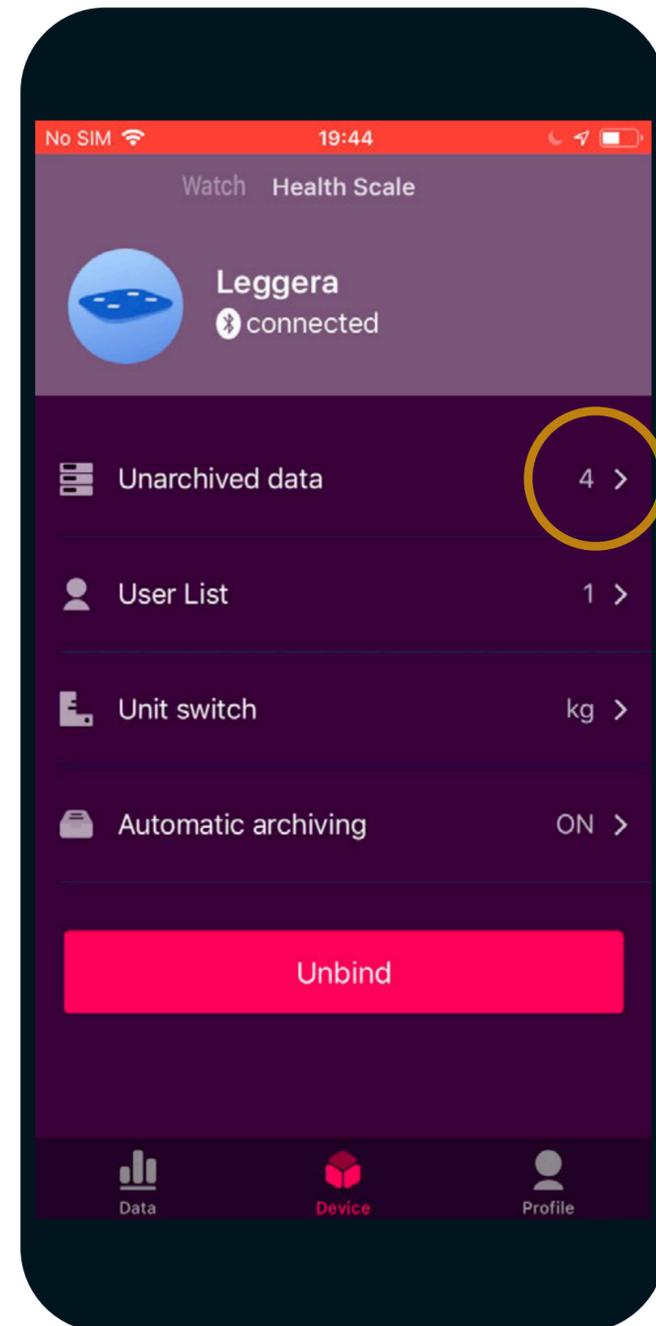



D. Tap on "Unarchived data" to see all the measurements saved within the Leggera memory.

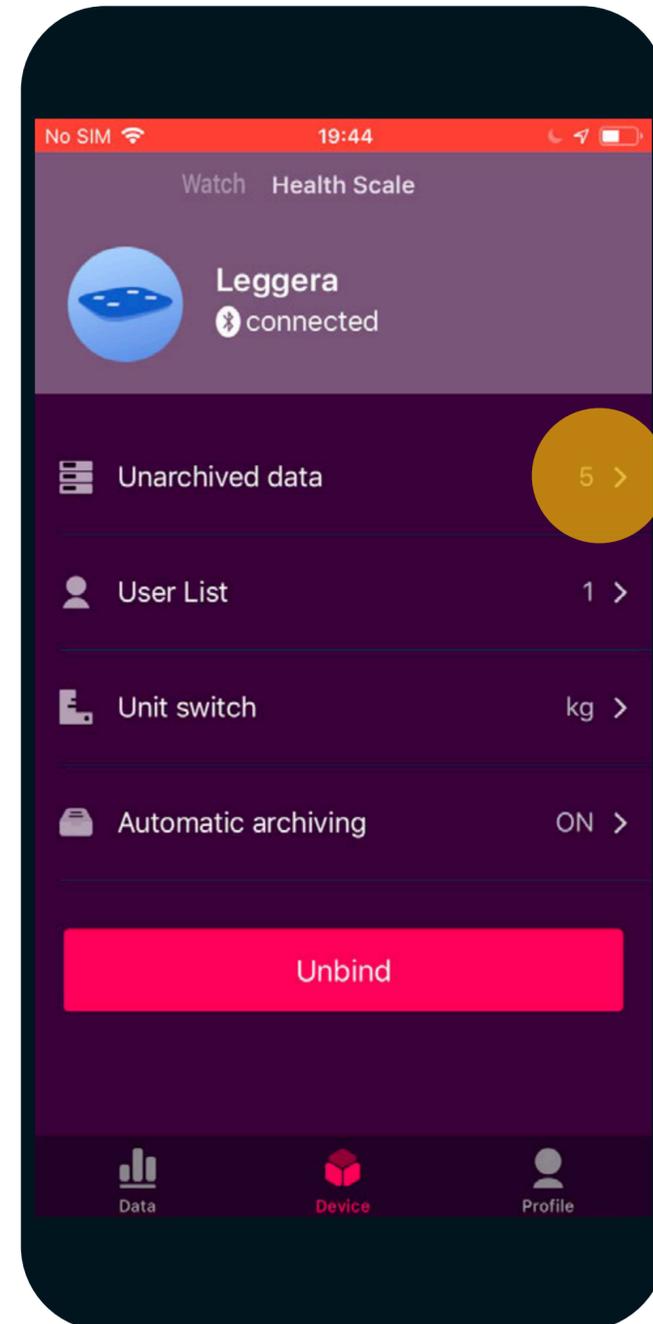


E. Tap on the last measure and assign it to yourself.

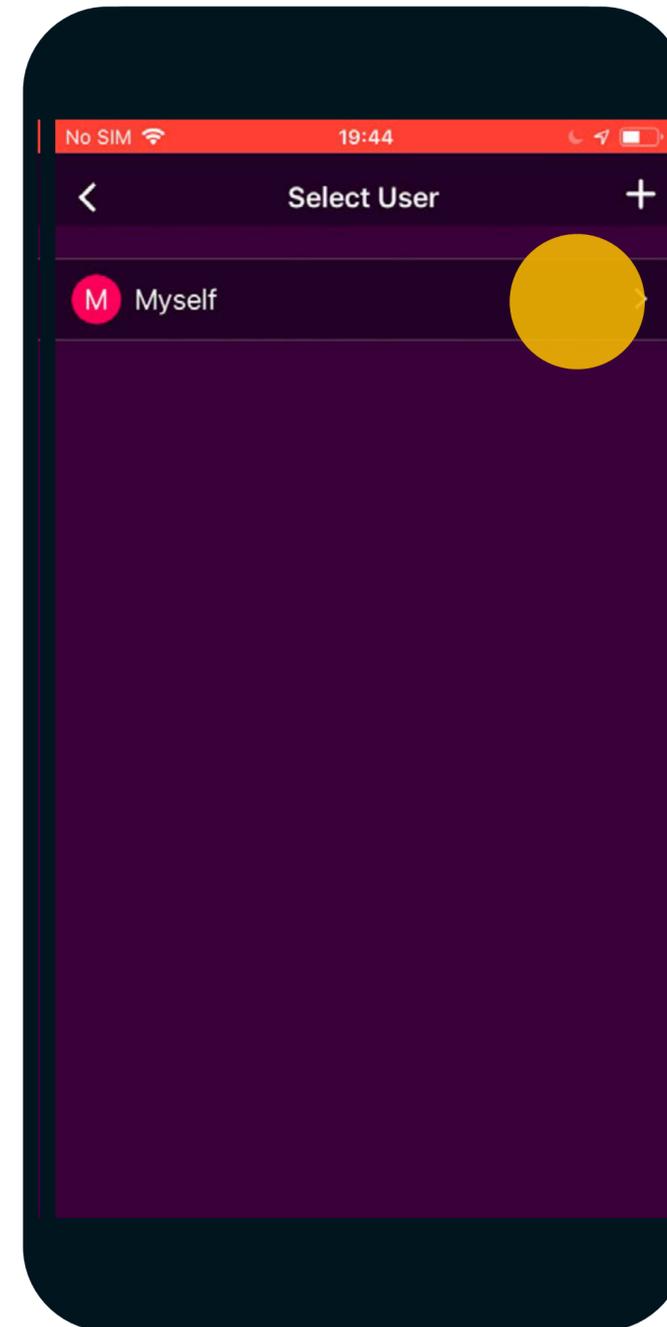




F.Leggera learns from you.
Simply step on Leggera and access the Device section.

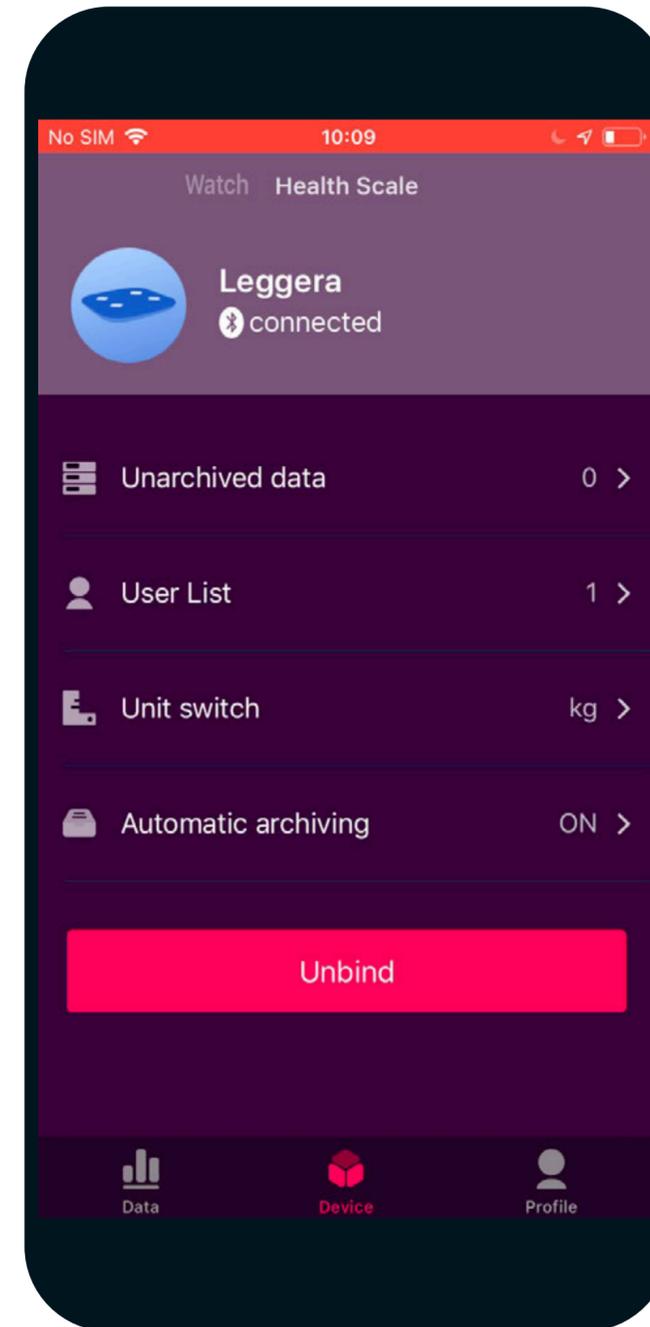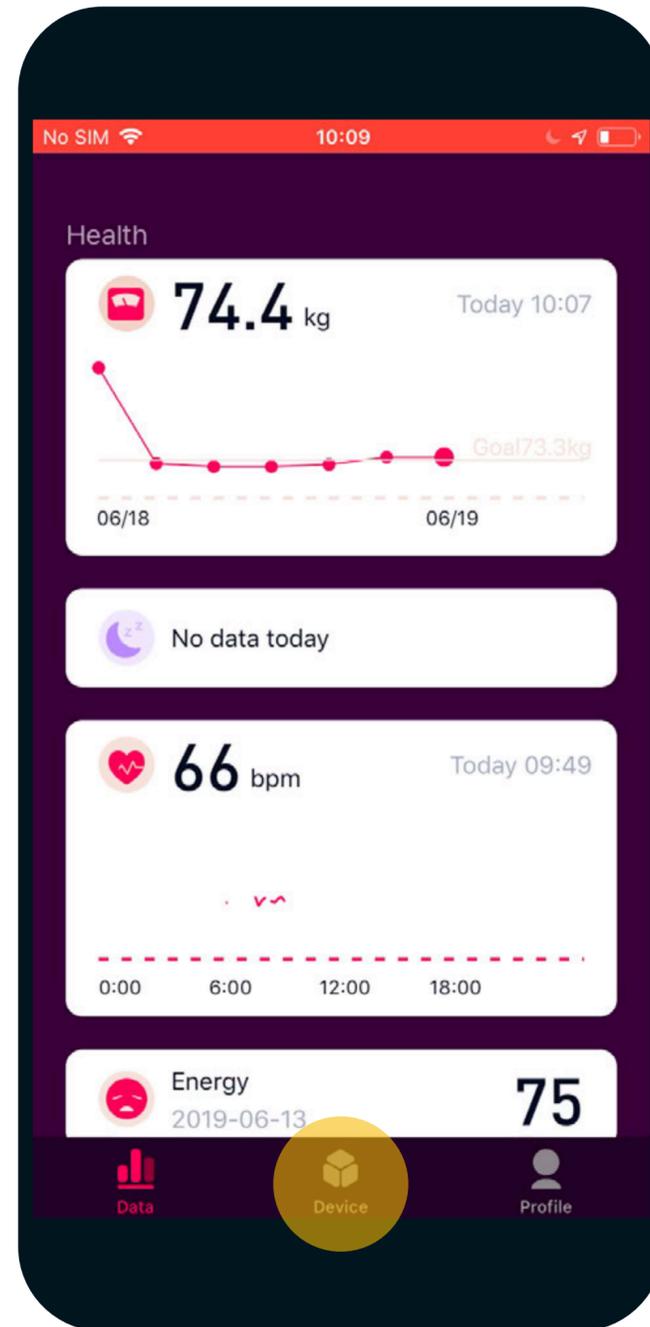

After a few measurements it will be able to identify the user and automatically assign the measure to the correct user.

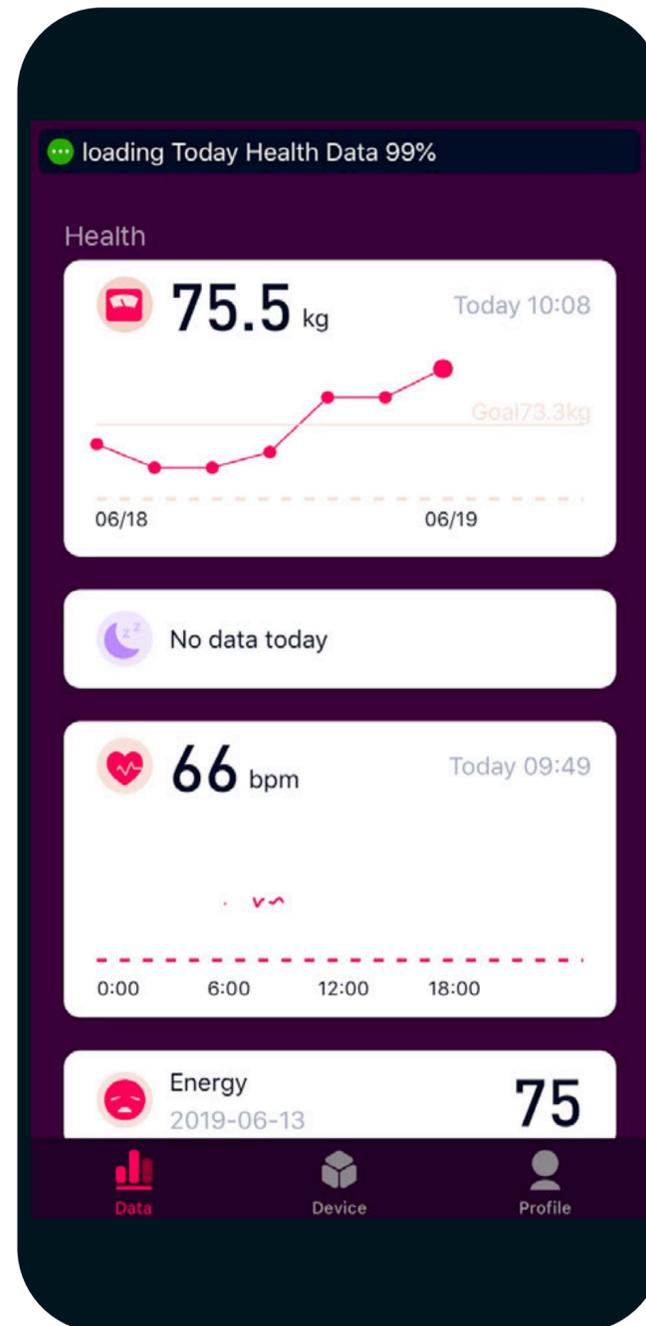

H. In the Home section, scroll to reach the health/weight section and open it.

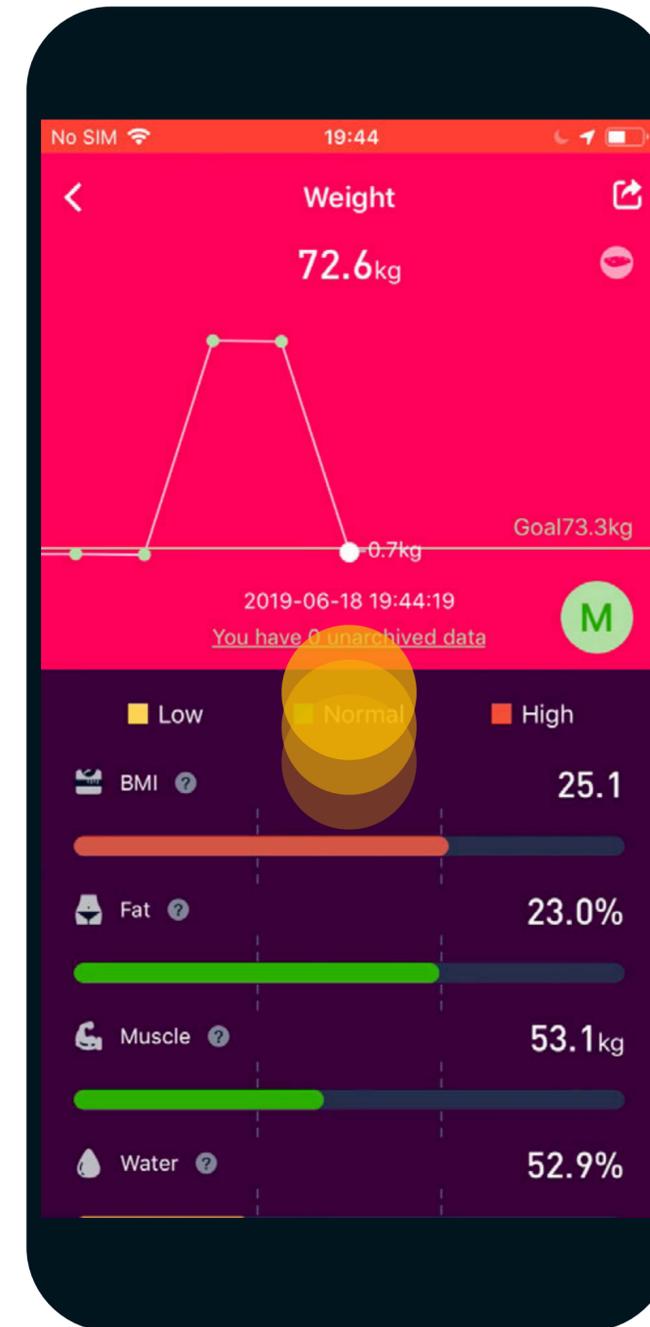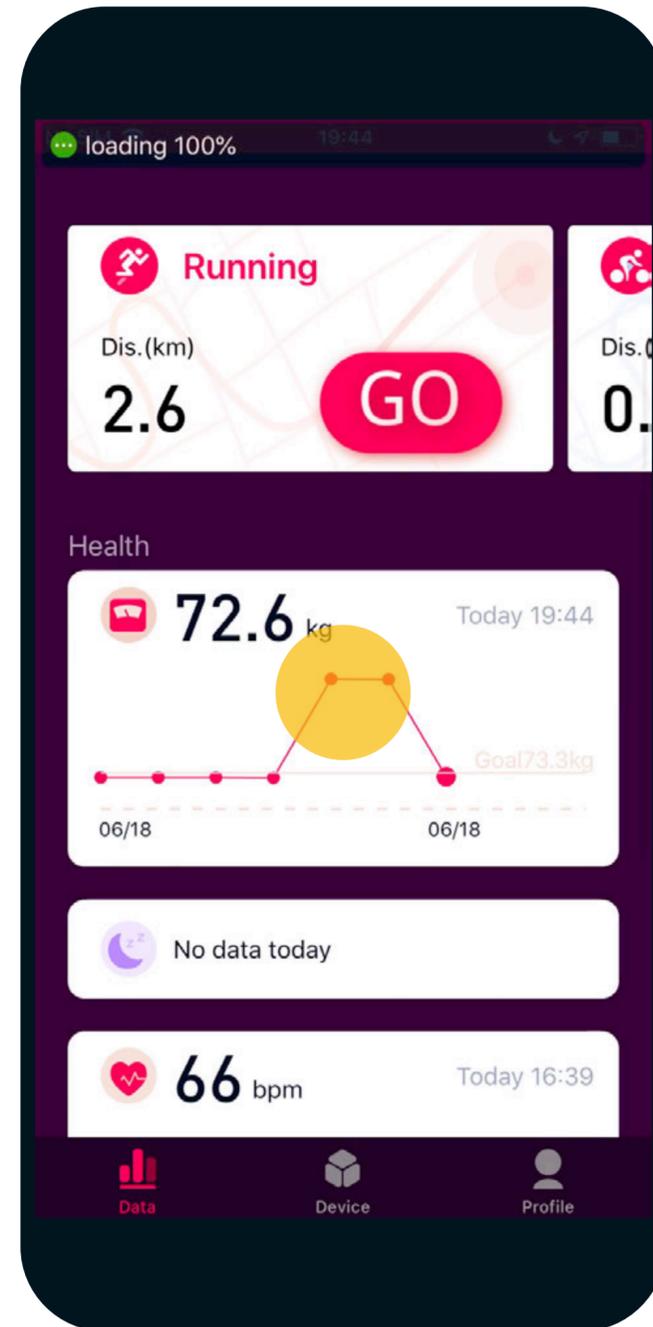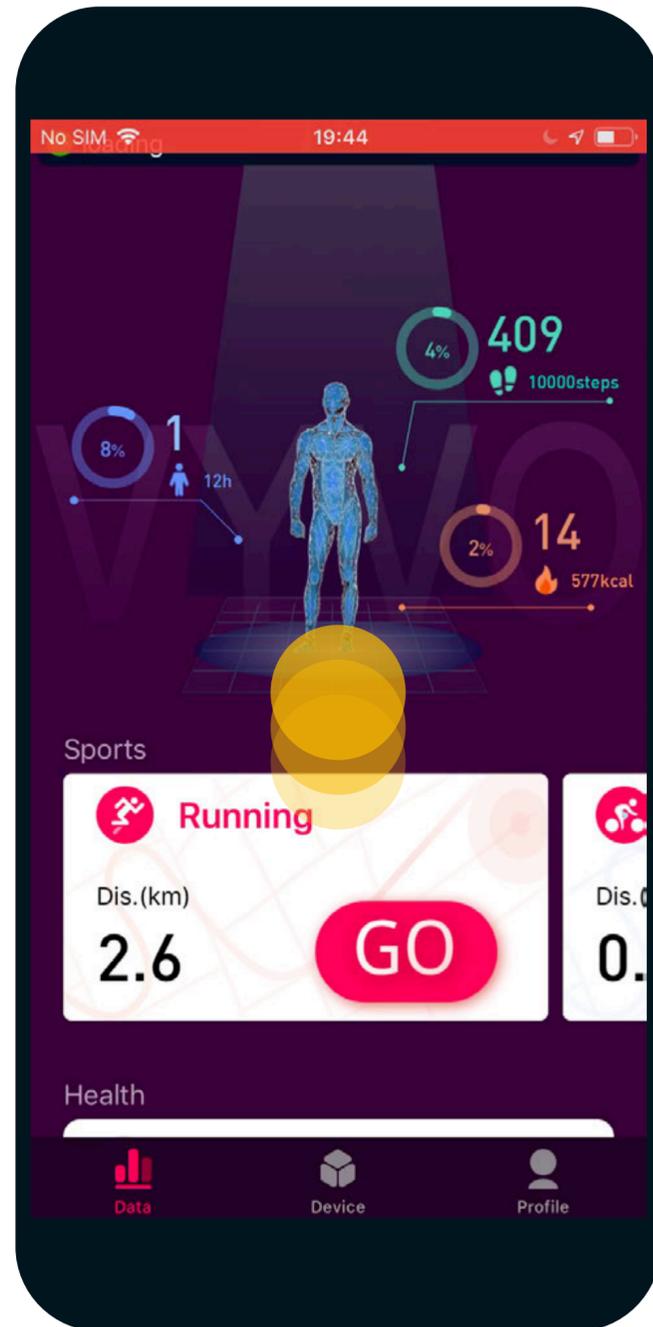

Here you can check all the details of your Leggera measurements.

vyvo™

MAKING LIFE WONDERFUL